

A woman with voluminous, curly dark hair is shown from the chest up, wearing a black, off-the-shoulder top. She is looking towards the camera with a slight smile. In her right hand, she holds a small, white, behind-the-ear hearing aid device. The background is plain white.

moxiTM

You can have it all
with a little Moxi.TM

MoxiTM style with YuuTM and NextTM

unitron[®]
hearing 


Look closely.

You'll love what you see.

Who couldn't use a little Moxi™? Moxi is the only discreet Canal Receiver Technology (CRT) hearing instrument able to meet the performance and style demands of clients with mild to severe hearing losses. Moxi's look is fashion-forward – but with such an ultra-small profile, you'll have to look closely to appreciate everything it has to offer:

- ▶ **Fits mild to severe hearing losses:** the most comprehensive fitting range within a single CRT product family with a choice of receivers
- ▶ **Control on demand:** onBoard™ control offers a choice of volume control or program button
- ▶ **OptimumFit™:** the industry's first adjustable custom CRT option, providing maximum comfort and the flexibility of adjustable tube lengths
- ▶ **Sleek and stylish:** 37% smaller than a traditional BTE
- ▶ **Premium plus to essential level:** Moxi style is available in fully featured Yuu™, Next™ 16, Next 8 and Next 4 hearing instruments



The image shows the Moxi hearing aid system components. On the left is the xReceiver unit, which includes a connector tube and a receiver. In the center is the main Moxi hearing aid, featuring the 'moxi' logo and an 'onBoard control' button. On the right is a close-up of the battery door, with a small circular button highlighted by two concentric orange circles. Labels with orange lines point to the receiver unit, a choice of couplings, and the onBoard control button.

xReceiver unit

includes:
connector tube
and receiver

choice of
couplings

onBoard™
control


**So much technology.
So little to see.**

Get onBoard™ – take control.

Most wearers prefer to retain control over the loudness of their hearing instruments. An industry first for new, smaller BTEs, Moxi features onBoard™ control: a push-button programmed within U:fit™ fitting software to be either a volume control or a program button. Incorporated discreetly into the battery door, onBoard control gives your clients control that hasn't been available in a discreet CRT, and allows you to tailor Moxi to individual preferences.

No compromise: automatic telephone solution in every Moxi.

With Moxi, there's simply no compromise in style or performance. In addition to Moxi's many advanced performance and comfort features, every Moxi-styled product includes the automatic phone solution with easy-t and telecoil.

A close-up of the Smart Control remote, showing several buttons with icons: a green checkmark, a hand, a plus sign, a circular arrow, and a square with a diagonal line. The text 'smart control' is visible on the device.


**So much more than your
average remote control:**

Moxi is available with the optional
Smart Control, offering added
functionality and convenience.

Unparalleled fitting flexibility.

Moxi offers a choice of receivers within one CRT product family, to satisfy clients with hearing losses ranging from mild to severe. A wide array of coupling and venting options ensures optimal performance for every client. Moxi brings instant fit capability or custom options to both standard and power fittings.

Standard xReceiver

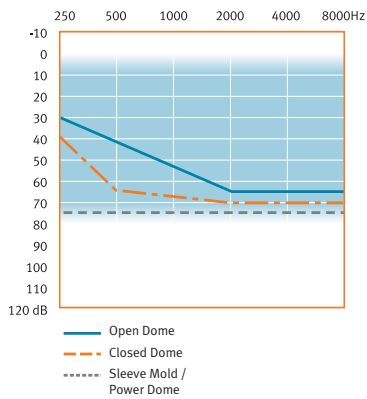

Open dome:
up to 65 dB HL

Choose for:

- ▶ Instant fit
- ▶ Open fit

Considerations:

- ▶ Occlusion free



Closed dome:
up to 70 dB HL

Choose for:

- ▶ Instant fit
- ▶ Semi-occluded fit

Considerations:

- ▶ Extended fitting range over open dome
- ▶ Some occlusion



Power dome:
up to 75 dB HL

Choose for:

- ▶ Instant fit
- ▶ Occluded fit

Considerations:

- ▶ Extended low-frequency response
- ▶ Occlusion



Sleeve mold:
up to 75 dB HL


Choose for:

- ▶ Custom fit
- ▶ Choice of venting

Considerations:

- ▶ Comfort
- ▶ Possibly occlusion, depending on venting
- ▶ Impression required

Power xReceiver

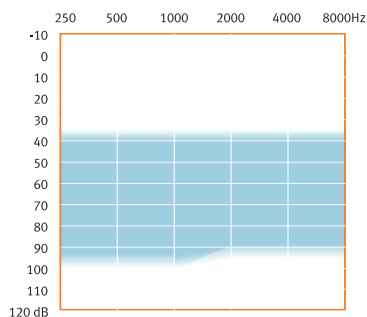

Power dome:
up to 95 dB HL

Choose for:

- ▶ Instant fit



Power shell:
up to 95 dB HL


Choose for:

- ▶ Custom fit
- ▶ Comfort
- ▶ Choice of venting

Considerations:

- ▶ Impression required

OptimumFit™

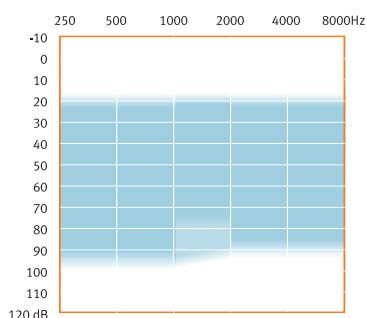


Custom shell:
up to 95 dB HL

Choose for:

- ▶ Custom fit
- ▶ Comfort
- ▶ Occlusion free
- ▶ Choice of venting

Considerations:

- ▶ Impression required
- ▶ Adjustable tube for optimized fit

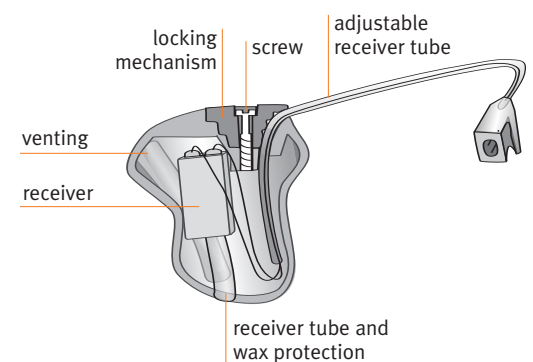

Moxi™ at-a-glance

Styled for:



premium plus

essential

Performance profile

Range of loss	Mild to severe	Mild to severe	Mild to severe	Mild to severe
Channels	20	16	8	4
Bands	20	16	8	8
Processing types	WDRC	WDRC Linear Limiting	WDRC Linear Limiting	WDRC Linear Limiting

Features

Automatic programs	AutoPro4™ (4 sound destinations)	AutoPro3™ (3 sound destinations)	AutoPro2™ (2 sound destinations)	
Manual listening programs	Up to 3	Up to 3	Up to 3	Up to 4
OnBoard™ control (VC or program control)	✓	✓	✓	✓
Directional microphone system	Multiband adaptive directionality	Adaptive directionality	Adaptive directionality	Fixed directional
Comfort-Clarity Balance	✓	✓		
AntiShock™	✓	✓	✓	✓
Speech enhancement LD	✓	✓	✓	✓
Phase canceller	✓	✓	✓	✓
Advanced Data logging	✓	✓	✓	✓
Self learning	✓			
Smart Control (optional)	✓	✓	✓	✓
LearnNow™	with Smart Control			
MyMusic™	✓	✓	✓	
Noise reduction	✓	✓	✓	✓
Wind noise manager	✓	✓	✓	✓
Easy-t	✓	✓	✓	✓
Telecoil	✓	✓	✓	✓
Program beep indicator	✓	✓	✓	✓
Low battery warning	✓	✓	✓	✓
Multiple choices of instant fit	✓	✓	✓	✓
Insitu Audiogram	✓	✓	✓	✓

Four products. Four price points.
One CRT product family.





Creating innovative hearing solutions that make everyday listening easier for people with hearing loss.