

Phonak Naída™ V

Product information



Phonak Naída V provides enhanced hearing performance for people with severe to profound hearing loss. With 4 performance levels (90, 70, 50, 30) and 3 models (RIC, SP, and UP), Naída V utilizes groundbreaking, dedicated features such as SoundRecover2 and BroadbandBooster for overall increased audibility and speech understanding. With more power than ever before, the fitting ranges of Naída V-RIC and Naída V-SP have been extended towards more severe losses allowing more people to benefit from a smaller device.

Naída V combined with Roger is the solution to effectively overcome listening challenges in difficult listening situations.




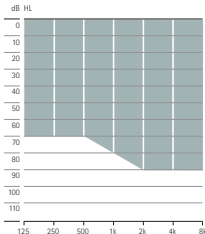
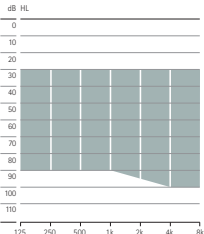
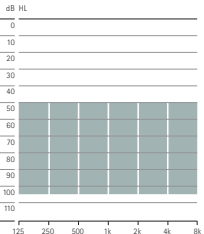


Call
Absolute Hearing Solutions
For Hearing Aid Deals
614-452-4280

www.absolutehearingsolutions.com

For more information visit www.phonakpro-us.com

PHONAK
life is on

Product description

							
	Naída V-RIC (V90/V70/V50/V30)			Naída V-SP (V90/V70/V50/V30)		Naída V-UP (V90/V70/V50/V30)	
Battery size	13			13		675	
Push button (Program/Volume)	•			•		•	
Volume control	•			•		•	
Telecoil	•			•		•	
EasyPhone	•			•		•	
Direct Audio Input	•			•		•	
Phonak CROS II compatibility	•			•		•	
IP rating	IP68 ¹			IP68 ¹		IP68 ¹	
Dimensions (L x W x D)	31.6 x 14.6 x 8.4 mm (1.24 x 0.57 x 0.33")			38 x 17.9 x 8.6 mm (1.50 x 0.70 x 0.34")		42.5 x 19.2 x 8.7 mm (1.68 x 0.76 x 0.34")	
Weight	2.2 g (0.08 oz)			4.6 g (0.16 oz)		5.4 g (0.19 oz)	
	xS	xP	xUP	HE11 680	HE11	HE11 680	HE11
Maximum Power Output (dB SPL)							
2 cc coupler	111	124	130	133	139	135	141
Ear simulator	122	132	136	139	143	140	145
Maximum gain (dB)							
2 cc coupler	46	57	66	74	81	77	82
Ear simulator	56	65	72	77	83	81	85
Frequency range – Ear simulator (Hz)	<100 - 9200			<100 - 5000		<100 - 5100	
Fitting range							
	 Power SlimTube II			 Hook/mini hook			
	Mild to moderately-severe hearing loss, all audiometric configurations			Mild to severe hearing loss, all audiometric configurations.		Moderate to profound hearing loss, all audiometric configurations.	

¹ IP68 indicates that the hearing aid is water resistant and dust tight. It survived continuous immersion in 1m of water for 60 minutes and 8 hours in a dust chamber as per the IEC60529 standard. No traces of dust were evident within the housing.

Housing colors

Hair & Skin



P1
Sand Beige



P2
Amber Beige



P3
Sandalwood



P4
Chestnut



P5
Champagne



P6
Silver Gray



P7
Graphite Gray



P8
Velvet Black

Fashion



P9
Ruby




Q1
Petrol

Traditional



O1
Beige













Acoustic couplings

				Naída V-RIC	Naída V-SP	Naída V-UP
Sound outlet	Hook standard	Damped	Transparent		HE11 680	HE11 680
		Undamped	Transparent		HE11	HE11
	Power SlimTube II	Left and right, lengths (00, 0, 1, 2, 3)			•	•
Receiver Length		Left and right, lengths (0, 1, 2, 3)		•		
Earpiece for xS and xP receiver		Dome	Open Dome (S, M, L)	•		
			Closed Dome (S, M, L)	•		
			Power Dome (S, M, L)	•		
SlimTip for xS and xP receiver		Hard material (only option for xP receiver)	Hollow shell	•		
		Soft material	Solid shell	•		
cShell for xS, xP and xUP receiver		Hard material (only option for xUP receiver)		•		
		Soft material		•		

Fitting








Software	Phonak Target™ 4.3 or higher
Interfaces	iCube II, NOAHLINK, HI-PRO, HI-PRO2

Roger Portfolio

Roger receivers		Naída V-RIC (V90/V70/V50/V30)	Naída V-SP (V90/V70/V50/V30)	Naída V-UP (V90/V70/V50/V30)
	Roger 18	•	•	
	Roger 19			•
	Roger X / AS18	•	•	
	Roger X / AS19			•
	Roger X / Phonak ComPilot II	•	•	•
	Roger MyLink	•	•	•
Roger microphones		Naída V-RIC (V90/V70/V50/V30)	Naída V-SP (V90/V70/V50/V30)	Naída V-UP (V90/V70/V50/V30)
	Roger Pen	•	•	•
	Roger EasyPen	•	•	•
	Roger Clip-On Mic	•	•	•
	Roger Table Mic	•	•	•

For information on Roger microphones please visit www.phonakpro-us.com/roger

Wireless accessories

		Naída V-RIC (V90/V70/V50/V30)	Naída V-SP (V90/V70/V50/V30)	Naída V-UP (V90/V70/V50/V30)
	Phonak DECT II	•	•	•
	Phonak PilotOne II	•	•	•
	Phonak EasyCall II	•	•	•
	Phonak ComPilot II	•	•	•
	Phonak ComPilot Air II	•	•	•
	Phonak TVLink II	• ¹	• ¹	• ¹
	Phonak RemoteMic	• ¹	• ¹	• ¹

¹ Requires Phonak ComPilot II or ComPilot Air II

Dedicated severe to profound features and benefits

Speech Understanding

SoundRecover2: The next generation of SoundRecover further improves audibility and utilizes an adaptive frequency lowering algorithm. The sound quality of low and mid frequencies are maintained while high-frequency sounds are made audible. This creates an optimal listening experience for more users.

Roger and directional setting: For the first time in any device, adults now have access to the proven benefits of directional microphones throughout their day when using Roger. In Phonak Target, the "Roger + Mic" program offers a choice of three microphone modes: omnidirectional, Real Ear Sound or an automatically activated fixed directional.

Audibility

BroadbandBooster: A new way of digitizing, sampling and filtering signals, boosts gain and output for loud broadband input signals, like speech. Mid and loud incoming broadband signals are more audible.

New mechanical and acoustical design: The new **xUP receiver** for Naída V-RIC offers additional gain and output with up to 6 dB in the low frequencies without increased battery consumption.

Robustness

With a new **glass fiber reinforced material**, Naída V is now 60% more robust than its predecessor which gives the client optimal confidence.

Water and dust resistant with an **IP68 rating**. IP68 indicates that the hearing aid is water resistant and dust tight. It survived continuous immersion in 1m of water for 60 minutes and 8 hours in a dust chamber as per the IEC60529 standard. No traces of dust were evident within the housing.

Size

Extended fitting ranges: the fitting ranges of Naída V-RIC and Naída V-SP have been extended towards more severe losses, allowing more people to benefit from smaller devices.

25% thinner than their predecessors, Naída V-SP and Naída V-UP are more comfortable to wear and less visible.

Simplicity

The new **RogerReady** feature in Naída V automatically activates the "Roger + Mic" program without the need for any extra programming of the hearing aid, so a Roger receiver can be added at any time, hassle-free.

Performance levels

	Premium (V90)	Advanced (V70)	Standard (V50)	Essential (V30)
AutoSense OS	Premium	Advanced	Standard	Essential
Calm Situation	•	•	•	•
Speech in Noise	•	•	•	•
Comfort in Noise	•	•	•	
Music	•	•		
Speech in Loud Noise	•			
Speech in Car	•			
Comfort in Echo	•			
Additional programs				
Max. additional programs	5	4	3	2
Speech in Wind	•			
Comfort in Echo	•			
Speech in Loud Noise	•	•		
Speech in 360°	•	•		
Speech in Noise	•	•	•	•
Calm Situation	•	•	•	•
Comfort in Noise	•	•	•	•
Music	•	•	•	•
Acoustic phone	•	•	•	•
Custom program	•	•	•	•
Streaming programs				
Max. streaming programs	4	4	3	3
Bluetooth audio + mic	•	•	•	•
Bluetooth phone / DECT + mic	•	•	•	•
Roger / RemoteMic	•	•	•	•
Audio jack	•	•	•	•

	Premium (V90)	Advanced (V70)	Standard (V50)	Essential (V30)
RogerReady	•	•	•	
Roger and directional setting	•	•	•	
UltraZoom	Premium	Advanced	Standard	Essential
SNR-Boost	•	•	•	
FlexControl	•	•	•	
FlexVolume	•	•	•	
DuoPhone	•	•	•	
SoundRecover2	•	•	•	•
BroadbandBooster	•	•	•	•
User Preference Tuning	•	•	•	•
Real Ear Sound	•	•	•	
Fine Tuning Channels	20	16	12	8
WhistleBlock	•	•	•	•
NoiseBlock	•	•	•	•
WindBlock	•	•		
EchoBlock	•			
SoundRelax	•	•		
Tinnitus Balance noise generator	•	•	•	•
QuickSync	•	•	•	•
auto Acclimatization	•	•	•	•

Life is on

We are sensitive to the needs of everyone who depends on our knowledge, ideas and care. And by creatively challenging the limits of technology, we develop innovations that help people hear, understand and experience more of life's rich soundscapes.

Interact freely. Communicate with confidence.
Live without limit. Life is on.

www.phonakpro-us.com



Absolute Hearing Solutions LLC
absolutehearing@att.net

750 Cross Pointe Road Suite F
Gahanna, Ohio 43230

Phone: 614-452-4280
Toll Free: 888-803-2159
Fax: 614-577-0481

www.absolutehearingsolutions.com